happybrush

SCHALL-ZAHNBÜRSTE

SONIC TOOTHBRUSH | BROSE À DENTS SONIQUE
SPAZZOLINO DA DENTI SONICO | SONISCHE TANDENBORSTEL

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DEUTSCH

3

ENGLISCH

15

FRANZÖSISCH

27

ITALIENISCH

39

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51

HALLO, ICH BIN HAPPYBRUSH!

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Hotline: 089 24412635
E-Mail: support@happybrush.de
Web: www.happybrush.de

Importeur: DS Direct GmbH | Am Heisterbusch 1, 19258 Golln

**PRODUKT**

happybrush
Schall-Zahnbürste HS03

KAUFDATUM

ZHB M DSDC V1 0717

TABLE OF CONTENTS

GERMAN 3

ENGLISH 15

FRENCH 27


ITALIAN 39

DUTCH 51

HI THERE, I'M HAPPYBRUSH!

For all you early birds, world travellers and night owls: I'm your new toothbrush and I'm here for you! **#gethappybrush**


OVERVIEW

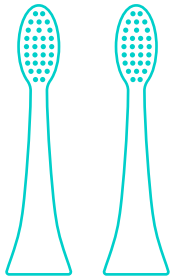2x Brush head

1x Charging station

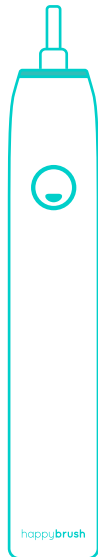1x Electric toothbrush

SAFETY INSTRUCTIONS

EXPLANATION OF SYMBOLS

Safety instructions: Read through these carefully and follow them to prevent personal injury and damage to property.

Protection class II and III

- This device may be used by children and by people with reduced physical, sensory or mental abilities or a lack of experience and/or knowledge if they are supervised or have received instruction on how to use the device safely and have understood the dangers resulting from failure to comply with the relevant safety precautions.
- Cleaning and user maintenance may not be performed by children, unless they are supervised.
- Children should not be allowed to play with the device.
- The charging station must always be disconnected from the mains power when

left unattended and before cleaning it.

- Always check the charging station for damage before putting it into operation. Only use it if the charging station, mains cable and mains plug have no signs of damage!
- If the electric toothbrush, the charging station, the mains cable or the mains plug of the device is damaged, it must be replaced by the manufacturer, customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid any hazards.

Danger of Electric Shock from Water!

- Never immerse the charging station, mains cable or mains plug in water or other liquids and ensure that these components cannot fall into water or become wet.
- Use the charging station only in closed rooms.
- Should the charging station fall into water, cut off the power supply immediately. Never attempt to pull the charging station out of the water while it is still connected to the mains power!
- Never touch the charging station, mains cable or mains plug with wet hands when these components are connected to the mains power.

Risk of Damage to Material and Property

- Position the mains cable in such a way that it does not become squashed or bent and does not come into contact with hot surfaces.
- Remove the mains plug from the plug socket if a fault occurs during charging or before a thunderstorm.
- When you pull it out of the plug socket, always pull on the mains plug and never the mains cable!

Warning of fire!

- Do not use the device in rooms containing easily ignitable dust or poisonous and explosive fumes.
- Do not cover the charging station during charging.

The Proper Handling of the Battery

- If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you touch battery acid, rinse the affected area with abundant fresh water immediately and contact a doctor as soon as possible.
- The battery must not be taken apart, thrown into a fire, immersed in liquids or short-circuited.

MY FEATURES

EFFECTIVE VIBRATIONS

My sonic technology ensures highly effective vibrations and gentle cleaning.

HIGH-QUALITY BRUSH INDICATOR

As soon as my smile fades, it's time for a new brush head. So, as soon as I'm not looking happy, change me.

HIGH-PERFORMING POWER BATTERY

My power battery will be with you during your adventures up to 3 weeks.

CUSTOM CLEANING MODES

Bye-bye brushing by hand model! It was nice, but now I have four new options: Normal, Sensitive, Polish, Massage.

PRACTICAL "LOCK" FUNCTION

You can simply "lock" me when you travel (5 second power button).

INTEGRATED 2 MINUTE TIMER

I have a practical built-in timer that lets you know when the 2 minutes are over.

PREPARATION AND INSTRUCTIONS

CHARGING

1. Plug your charging station into the power mains.
2. Pop your happybrush onto the charging station. A complete charge of the battery takes approx. 16 hours. When the LED lights up red, your toothbrush is charging. When charging is complete, the LED turns a steady green. It is normal for the handle to become warm during charging.

Important: Charging is done inductively and can be affected by metallic objects near the charging station. Keep at least 10 cm between the charging station and any metallic objects.

BRUSH HEADS

Push the brush head onto the metal tip of the handle until it sits firmly.


BRUSHING INSTRUCTIONS

Moisten my bristles with a bit of water and then apply toothpaste - happybrush toothpaste works best. Now you can start brushing!

Hold the brush head at a 45-degree angle toward your gum line. When brushing your teeth, move the brush back and forth gently, with a sweeping movement over the front, back and top (chewing surface) of your teeth.

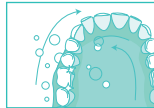
Start with the top back teeth and move to the front. Then do your bottom teeth, starting again at the back.

At the end of 2 minutes, the regular cleaning cycle is finished and you can polish areas with discolouration.

MODE OF OPERATION

CHARGING

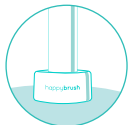While I am being charged, my LED lights up red. As soon as I am fully charged, it turns green.

READY, SET...

Use the power button to start me and change operating modes simply by pressing the button again (within 10 seconds).

GO!

I vibrate every 30 seconds so that you treat your 4 quadrants equally.

ALREADY DONE

After 2 minutes, I make two brief vibrations. You can keep on brushing if you want.

ON THE ROAD

If you press the power button for 5 seconds, I'll "lock" for travel.

CLEANING

HANDLE

1. Take the brush head off and rinse the metal tip with warm water. Please do not try to clean around the seal of the metal tip with a sharp object because it could damage the seal.
2. Now, simply wipe the handle with a dry cloth.

BRUSH HEAD

1. Rinse the brush head and bristles under running water after every use. Then just shake excess water off the brush head.
2. Remove the brush head from the toothbrush at least once a week and clean the connection between the handle and the brush head with warm water.
3. Use a new brush head every 2-3 months for best results. The indicator will tell you when: Simply change it when the smile fades.

CHARGING STATION

1. Remove the plug from the mains socket before cleaning the charging station.
2. Then wipe the charging station with a slightly damp cloth.

Important: The handle, the brush heads, and the charging station are not to be cleaned in the dishwasher. Do not submerge the charging station in water. Do not use abrasive cleaning agents. Also do not try to clean the charging station with any kind of sharp metallic object or nylon bristles.

DISPOSAL

RECYCLING

The environment is very important to us. When you no longer use your happybrush, do not throw it in with your normal household waste.

Dispose of your happybrush at a recycling centre for used electrical and electronic devices. More information is available from your local authorities.

To protect the environment, non-rechargeable and rechargeable batteries must not be disposed of with normal household waste, but must be taken to suitable collection points. Also note the relevant statutory regulations governing the disposal of batteries.

TECHNICAL SPECIFICATIONS

Power	3.0 W
Voltage	100-240 V - 50/60 Hz
Battery	3.7 V 800 mAh Li-Ion
Protection class	Class III ⚡ (Handle) Class II ⚡ (Charging station)
Protection level	IPX5 (Handle) IPX7 (Charging station)
Vibration	Normal: 30.000 per minute Sensitive: 26.000 per minute Polishing: 30.000 per minute Massage: 26.000-37.500 per minute
Type	HS03
Timer	4-quadrants timer (every 30 seconds)
Operation time	approx. 84 minutes
Charging time	Charge for 24 hours the first time, thereafter 16 hours

CE The CE label on this product is our way of confirming that it complies with the valid, product-specific European guidelines.

The voluntary certification "GS" stands for "Geprüfte Sicherheit" (certified safety) following the requirements of the product safety law. It certifies products that pose no risk to the user when used in compliance with operating instructions. Compliance is reviewed by independent testing companies such as TÜV SÜD.

WARRANTY

Your happybrush comes with a 1 year warranty as of the date of purchase. If defects of material or function should occur during the 1 year warranty period, we will repair or replace your happybrush.

The original proof of purchase (original receipt or payment confirmation) is to be sent to us as the warranty and proof of ownership. The receipt must show the date of purchase, your name and the full product type description of the device. The warranty applies only to the original purchaser. Should the device be sold to a second buyer, the warranty is no longer valid.

All claims from this warranty (which includes both the primary claim as well as any secondary claims after a repair or a replacement) are subject to a 12 month limitation period. The limitation period for the primary warranty starts from the moment the defect is known or gross negligent ignorance of the defect. The limitation period for secondary claims begins with the delivery of the repaired or replaced product.

The product may be replaced with the same model or a comparably equipped model. The replacement or repair of the product does not justify a new warranty. The warranty does not apply to devices or device parts subject to normal wear and tear and thus viewed as wear parts. The warranty is void if the device has been damaged, not used properly or not maintained properly. Repairs may only be done by happybrush GmbH. To ensure proper utilization of the device, all of the instructions which are listed in the operating instructions must be complied with precisely. Using it and handling it in ways that are not recommended in the instructions or for which warnings exist are to be absolutely avoided.

Important: The warranty and warranty period for specific devices and materials can differ from those shown here. These warranty restrictions have no impact on the statutory claim to warranty. The following device parts do not fall under the guarantee: Brush heads; damage due to misuse, wrong use, negligence or changing the product; normal wear and tear such as chipped spots, scratches, abrasions, discolourations or fading.

SUPPORT

happybrush GmbH
Amalienstr. 71
80799 Munich

Hotline: +49 (0)89 24412635
Email: support@happybrush.de
Web: www.happybrush.de

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